

MindWare 8-Slot BioNex Hardware

User Reference Guide

Copyright 2011 © by MindWare Technologies LTD. All Rights Reserved.



Copyright 2008 © MindWare Technologies Ltd.

MindWare 8-Slot BioNex Hardware User Reference Guide



Internet Support

E-mail: sales@mindwaretech.com

Web Address: www.mindwaretech.com

Phone

Ph: (614) 626-4888

Fax: (614) 626-4915

MindWare Technologies LTD.

1020-F Taylor Station Rd. • Gahanna, Oh 43230

Copyright 2011 © by MindWare Technologies LTD. All Rights Reserved.

Safety

Warning Regarding Medical and Clinical Use of MindWare Technologies LTD Products.

MindWare Technologies LTD products are not designed with components and testing for a level of reliability suitable for use in treatment and diagnosis of humans or as critical components in any life support systems whose failure to perform can reasonably be expected to cause significant injury to a human. Applications of MindWare Technologies LTD products involving medical or clinical treatment can create a potential for death or bodily injury caused by product failure, or by errors on the part of the user or application designer. Any use or application of MindWare Technologies LTD products for or involving medical or clinical treatment must be performed by properly trained and qualified medical personnel, and all traditional medical safeguards, equipment, and procedures that are appropriate in the particular situation to prevent serious injury or death should always continue to be used when MindWare Technologies LTD products are being used. MindWare Technologies LTD products are NOT intended to be a substitute for any form of established process, procedure, or equipment used to monitor or safeguard human health and safety in medical or clinical treatment.

Table of Contents

MindWare 8-Slot BioNex Hardware	1
Table of Contents.....	4
Section 1: Technical Specifications	5
Section 2: Instructions For Use.....	10
Identification of Operating Controls	10
Warnings.....	10
Accessories and Consumables.....	11
Section 3: Input/Output Connections	13
Back Panel	13
Front Panel.....	17
Section 4: Cleaning Instructions.....	19
Section 5: Additional/Optional Parts and Services.....	20
Module Installation Instructions.....	22
Section 6: Test System and Configuration.....	23

Section 1: Technical Specifications

Note: All timings are specified relative to the analog input stream

*** Specifications tested with minimum test system and configuration. Performance may vary with system configuration and load.**

Range of environmental conditions for which the unit is designed:

- Indoor use
- Altitude to 2000 m
- Temperature for safe operation: 5°C to 40°C
- Maximum relative humidity: 80% for temperature up to 31°C decreasing linearly to 50% at 40°C
- Power supply voltage fluctuations not to exceed $\pm 10\%$ of the nominal voltage
- Over voltage category II
- Pollution degree 2

Power Input:

- AC Input: 100-240VAC, 50-60Hz, 1.5A
- Fuse: 2A 250V 5x20mm Slo-Blow

Bipolar Power Output¹:

- $\pm 5V$, 90mA
- $\pm 12V$, 25mA
- Protection: External Fuse – 0.5A 250V 5x20mm

Unipolar Power Output²:

- 5V
- 300mA

MindWare 8-Slot BioNex Hardware User Reference Guide

- Protection: Auto-Reset Fuse

Analog Inputs:

- 32 Channels
- Resolution: 16-Bit
- Maximum Sample Rate: 250 kS/s over all channels
- Timing Accuracy: 50ppm of sample rate
- Range, accuracy, and noise depends on amplifier module

Digital Inputs:

- 8 Synchronous Inputs / 16 Synchronous Inputs²
 - Maximum Sample Rate: Synchronous with analog sample rate
 - Minimum Pulse Width: 2/Sample Rate
 - Minimum Time Between Pulses: 2/Sample Rate
 - Timing Accuracy: +/- 1/Sample Rate
 - Input High: 2.2V min, 5.25V max
 - Input Low: 0V min, 0.8V max
 - Input High Current: 250 uA
 - Input Low Current: -10uA
 - Pull Down Resistor: 50kohm typ., 20kohm min
- 8 Asynchronous Inputs¹ and Rating Response Buttons
 - Minimum Pulse Width: 70 ms *
 - Minimum Time Between Pulses: 100 ms *
 - Timing Accuracy: +/- 70 ms *

- Input High: 2.2V min, 5.25V max
- Input Low: 0V min, 0.8V max
- Input High Current: 250uA
- Input Low Current: -10uA
- Pull Down Resistor: 50kohm typ., 20kohm min

Trigger:

- 1 Trigger Input

- Trigger Accuracy: <1ms
- Pull Up Resistor¹: 4.7kohm
- Input High: 2.2V min, 5.25V max
- Input Low: 0V min, 0.8V max
- Input High Current: 250uA
- Input Low Current: -10uA

- 1 Trigger Output

- Maximum Propagation Delay from Trigger In on CE Marked Units: 250ns
- Maximum Propagation Delay from Trigger In on non-CE Marked Units: 10ms
- Pulse Width²: 50ms
- Output High: 2.5V min, 5.V max
- Output Low: 0V min, 0.5V max
- Output High Current: 8mA
- Output Low Current: -8mA

MindWare 8-Slot BioNex Hardware User Reference Guide

Analog Output:

- 2 Digital to Analog converters
- Resolution: 16-bit
- Range: +/- 10V
- Output Impedance: 0.2ohm
- Output Drive: +/- 5mA

Video Capture Module:

- Timing Accuracy: +/- 10 frames *
- Mode: NTSC or PAL
- Video Input: Composite (BNC), 75ohm
- Audio Input: Line In Stereo (3.5mm jack)
- Resolution: 360x240, 720x480
- Formats: MPEG-1, MPEG-2, MPEG-4, MJPEG

Legacy Video Capture Module:

- Timing Accuracy: +/- 10 frames *
- Mode: NTSC or PAL
- Video Input: Composite (BNC) or S-Video (DIN), 75ohm
- Audio Input: Line In Stereo (2x RCA), or Microphone mono (3.5mm 3-pos jack), 10kOhm
- Resolution: 360x240, 720x480

- Formats: MPEG-1, MPEG-2, MPEG-4, MJPEG

¹Feature only available in non-CE marked units

²Feature only available in CE marked units

Section 2: Instructions For Use

Identification of Operating Controls

The operating controls for this hardware are located in the software.

Warnings

Review the following safety precautions to avoid injury and prevent damage to the BioNex hardware and any products connected to it. To avoid potential hazards, use this equipment only as specified.



This symbol indicates the equipment is to be serviced by trained personnel only. Failure to comply risks electric shock.



WARNING: This symbol on the equipment indicates that more detailed information will be found in the instruction manual.



WARNING: The ground pin on the Power Cable must remain intact. Failure to comply risks electric shock.

When positioning the BioNex 8-Slot chassis ensure the power switch is easily accessible.

If any equipment is used in a manner not specified by the manufacturer the protection provided by the equipment may be impaired.

Accessories and Consumables



8-Slot E-Prime Trigger Cable

Part #40-0027-00

Instructions:

- Turn power to the BioNex hardware off
- Connect the DB25 with the pigtail RCA connector to the Digital I/O 1 or Digital I/O 2 port
- Connect the RCA connector to the Trigger Input jack
- With a small slotted screwdriver, gently screw the DB25 connector into place
- Connect the 2nd DB25 connector to the parallel port of the computer with E-Prime installed



0.5A 250 V 5X20 mm Fast Acting Fuse

Part #55-0833-00

Instructions:

- Turn power to the BioNex hardware off
- Unplug the power cable from the BioNex hardware
- With a slotted screwdriver, turn the fuse holder and remove the fuse from the chassis
- Replace the fuse and re-insert the fuse into the fuse holder



MindWare 8-Slot BioNex Hardware User Reference Guide




2A T 250 V 5X20 mm Slow Blow Fuse

Part #55-0834-00

Instructions:

- Turn power to the BioNex hardware off
- Unplug the power cable from the BioNex hardware
- With needle nose pliers, gently remove the fuse drawer from the power entry module.
- Replace fuse
- Re-insert the fuse drawer into the power entry module

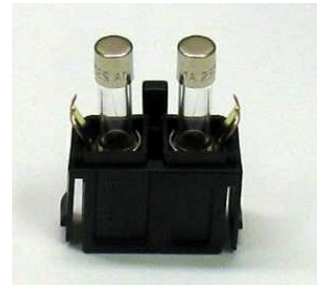

Section 3: Input/Output Connections

Back Panel



1

$\pm 5V$, 90mA

- Tip - +5V
- Ring - -5V
- Sleeve - Common

2

±12V, 25mA

- Tip - +12V
- Ring - -12V
- Sleeve - Common

3

TRIGGER INPUT – Event Trigger Input

- Tip – TTL Level In
- Sleeve - Ground

4

TRIGGER OUTPUT

- Tip – TTL Level Out
- Sleeve - Ground

5

DAC 1

- Digital to Analog Output – 1/8" Mono Audio Jack
 - Tip – DC Voltage out
 - Sleeve – Ground
 - 0-5V DC Range

6

DAC 2

- Digital to Analog Output – 1/8" Mono Audio Jack
 - Tip – DC Voltage out
 - Sleeve – Ground
 - 0-5V DC Range

7

FUSE

- Replace with 250V, 0.5A Fuse

8

SYNC

- Tip – Sync
- Sleeve - Ground

9

USB

- Standard USB Type-B connector

10

DIGITAL I/O 1

- Asynchronous Digital Input/Output
- Software Selectable
- DB25
 - Pin 1 - Pin 16 : Bit 0 – Bit 15
 - Pin 17 – Pin 25 : Ground

11

DIGITAL I/O 2

- Synchronous Digital Input/Output
- Software Selectable
- DB25
 - Pin 1 – Pin 16 : Bit 0 – Bit 15
 - Pin 17 – Pin 25 : Ground

12

POWER ENTRY MODULE

- AC Input: 100–240VAC, 50-60Hz, 1.5A
- Power Toggle Switch
- 2A T 250 V 5X20 mm Fast Acting Fuse

Front Panel



1

COMP

- Standard BNC Composite Video connection

2

S-VID

- Standard S-Video connection

3

AUDIO IN

- Left and Right RCA audio connections

4

MIC IN

- Standard 1/8" Audio connection

5

SUBJECT IN

- Amphenol PL-700 Connection
- Requires MindWare Subject Cable Dual 5 Position (Part #40-0020)

6

TRANSDUCER/RATING DIAL IN

- Amphenol PL-500 Connection
- Requires MindWare Transducers or MindWare Rating Dials to function properly

7

HIGH LEVEL IN

- Standard 1/8" Audio connection

Section 4: Cleaning Instructions

To remove scratches on front panels:

- With an eraser gently rub across the scratch
- Clean off any debris before using

Cleaning Hardware:


- Turn the power to the BioNex off
- Disconnect the power and all other cables connected to the device
- With a slightly damp cloth, wipe off hardware
 - o Avoid getting moisture in exposed connectors
- Allow to dry thoroughly before reconnecting power

Section 5: Additional/Optional Parts and Services

Rack Mounting

- Can be ordered with the BioNex
- Equipment must be returned to the manufacturer to have the Rack Mounting hardware installed.

Video Acquisition

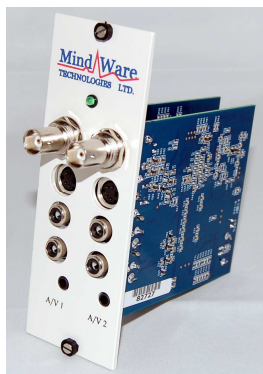

- Available in Single, Dual, and Quad modes
- Can be ordered with the BioNex hardware
- Equipment must be returned to the manufacturer to have the Video module(s) installed

Acquisition Modules



Impedance Cardiograph & GSC

- Optional with the BioNex hardware
- Equipment can be ordered separately from the BioNex hardware for installation in existing systems
 - See Module Installation below for installation instructions
- Subject harness sold separately (Part #40-0020)




3 Channel Bio-Amplifier & GSC

- Optional with the BioNex hardware
- Equipment can be ordered separately from the BioNex hardware for installation in existing systems
 - See Module Installation below for installation instructions
- Subject harness sold separately (Part #40-0020)



4 Channel Bio-Amplifier

- Optional with BioNex hardware
- Equipment can be ordered separately from the BioNex hardware for installation in existing systems
 - See Module Installation below for installation instructions
- Subject harness sold separately (Part #40-0020)

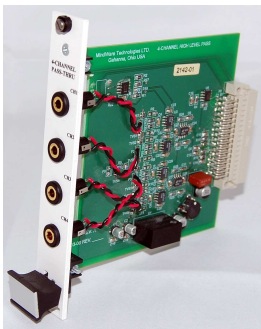

4 Channel High Level Pass Through


- Optional with BioNex hardware
- Equipment can be ordered separately from the BioNex hardware for installation in existing systems
 - See Module Installation below for installation instructions



4 Channel Rating Response

- Optional with BioNex hardware
- Equipment can be ordered separately from the BioNex hardware for installation in existing systems
 - See Module Installation below for installation instructions
-

MindWare 8-Slot BioNex Hardware User Reference Guide


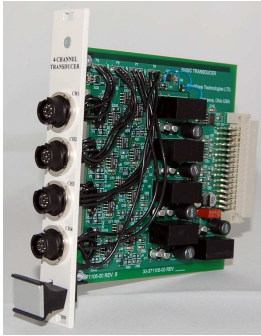


4 Channel Transducer

- Optional with BioNex hardware
- Equipment can be ordered separately from the BioNex hardware for installation in existing systems
 - See Module Installation below for installation instructions
- Transducers are sold separately

Module Installation Instructions

- Turn off the power to the BioNex hardware
- Unplug the power cable from the BioNex hardware
- Slide Module into empty slot
 - Make sure the card is in the top and bottom rails
- With a small slotted screwdriver gently screw the top and bottom screws into place

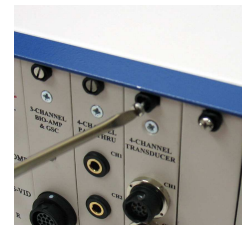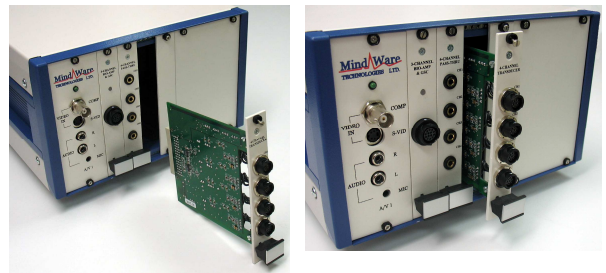

Section 6: Test System and Configuration

Minimum System Requirements:

- 1024x768 graphics resolution
- USB 2.0 with no USB hubs between machine and BioNex
- 4GB RAM for Windows XP
- 8GB RAM for Windows 7

Basic Test System:

- 2.8GHz Intel Core 2 Duo
- 4GB RAM
- Windows XP

BioNex 8-Slot Testing Configuration

- BioLab version 3.0.8
- 1000 Samples/sec
- 1/sec Update Rate
- 16 Analog Input channels enabled
- 8 Synchronous Digital Event channels enabled
- 8 Asynchronous Digital Event channels enabled
- 2 Channels of video at 720x480 resolution, MPEG-4 encoded

MindWare 8-Slot BioNex Hardware User Reference Guide

DECLARATION OF CONFORMITY

Application of Council Directive: Low Voltage Equipment Directive (2006/95/EC)

Standards to which Conformity is declared: EN 61010-1:2001

Application of Council Directive: Electro Magnetic Compatibility (2004/108/EC)

Standards to which Conformity is declared: EN 61326-1:2006

Manufacturer Name: MindWare Technologies Ltd.

Manufacturer Address: 1020-F Taylor Station Rd.

Gahanna, OH 43230

USA

Type of Equipment: Behavioral Research Monitor

Model Number: 50-3711-08

Serial Number(s) All First year of manufacture 2009

I, the undersigned, hereby declare that the equipment specified above conforms to the above Directive and Standards, and this Declaration is supported by a Technical File located at the Factory. Conformity Assessment is consistent with the requirements of the Directives.

Place - Gahanna, Ohio, USA

(Signature)

9-1-09
(Date)

Martin G. Gillman
(Printed Name)

V.P. Engineering
(Position)



1020-F Taylor Station Road

Gahanna, OH 43230

(614) 626- 4888

www.mindwaretech.com

Copyright 2011 © MindWare Technologies LTD.